

Actim[®] PROM

A highly specific test for detecting premature rupture of fetal membranes.

- Actim[®] PROM identifies levels of IGFBP-1 in vaginal swap samples to detect premature rupture of fetal membranes.
- Test results are available at the patient's bedside in just 5 minutes.
- Blood and other common contaminants do not interfere with the test results.

**No interference
with blood
or semen, urine,
vaginal medications,
lubricants, bathing
products or
infections.**



A pregnant woman with dark hair and a warm smile is looking towards a healthcare provider whose back is partially visible on the right. The woman is wearing a white tank top under a grey cardigan. The background is softly blurred, showing what appears to be a clinical or office setting.

Actim PROM is a **one-step** dipstick test that gives **results in just 5 minutes** with sampling completed in seconds – with or without a speculum.

Actim PROM can be used at **any gestational age**. Furthermore, blood, semen, urine, vaginal medications, lubricants or infections do not interfere with the test results.

About 20% of pregnant people with suspected PROM **have vaginal bleeding**. All these patients can be diagnosed with the help of Actim PROM.

Actim PROM is helping **millions of pregnant people** with suspected PROM around the globe.

Actim PROM test results **can be digitally interpreted** with Actim 1ngeni, taking data connectivity one step further in pregnancy monitoring.

A unique rapid test that detects IGFBP-1 levels in vaginal fluids

Actim PROM is based on highly specific and unique monoclonal antibodies that bind to the insulin-like growth factor binding protein-1 (IGFBP-1) present in amniotic fluid in high amounts.

The concentration of IGFBP-1 in amniotic fluid rises early in pregnancy and remains elevated until birth. Amniotic fluid is not normally found in the vagina, but when fetal membranes rupture, amniotic fluid leaks into the vagina and the IGFBP-1 concentration rises quickly. Actim PROM detects IGFBP-1 in vaginal swab samples and helps to identify premature rupture of membranes.

Thanks to its optimal detection limit, Actim PROM can identify even small amounts of amniotic fluid with a minimal risk of false positive results due to the very low levels of IGFBP-1 normally found in the vagina during pregnancy.

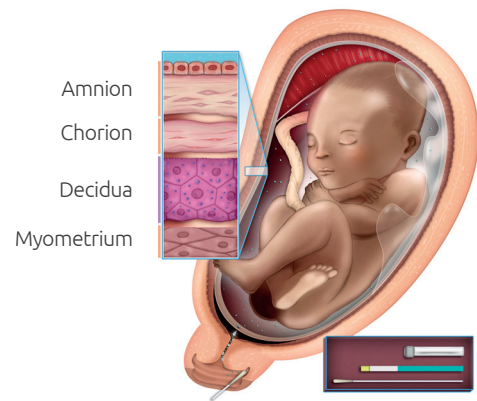


Figure 1. Actim PROM identifies premature rupture of membranes through a simple vaginal swab sample.

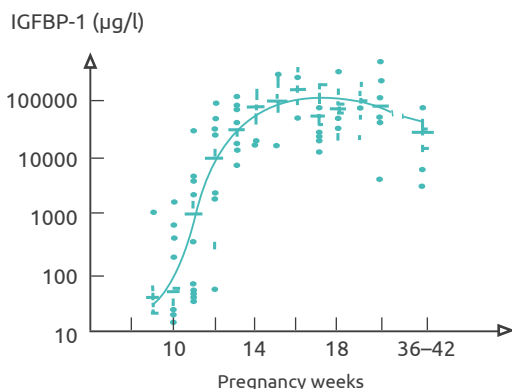


Figure 2. IGFBP-1 concentration in amniotic fluid rises quickly in early pregnancy and remains high until birth (Wathen et al. 1993).

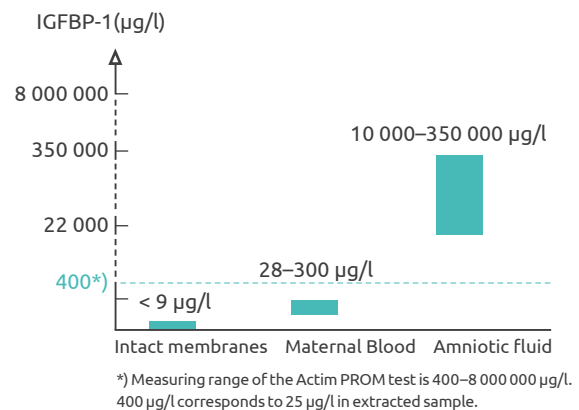


Figure 3. The detection range of Actim PROM covers all clinically relevant concentrations of IGFBP-1.

Complications of PROM

Premature rupture of membranes (PROM) occurs when the fetal membranes break before the onset of labor. Once the membranes break, both the birth-giver and the child are at high risk of infection and other complications.

Premature rupture of membranes causes complications in 2% to 20% of pregnancies. It can occur at any gestational age, and it eventually leads to the start of labor and delivery. PROM is a leading cause of approximately one third of all preterm births, and it is associated with one fifth of perinatal deaths.

A reliable test for all patients suspected with PROM

IGFBP-1-testing is mentioned in several national guidelines as a recommendation to detect PROM when the rupture of membranes remains uncertain after a full evaluation.*

Actim PROM is optimized to be **so sensitive that it detects even the smallest ruptures** that are clinically invisible (even less than 1 µl of amniotic fluid). These tiny ruptures cannot be detected with traditional methods, but are clinically relevant as they can induce delivery, cause infections, and threaten the health of both the birth-giver and the fetus.

Thanks to Actim PROM's **specificity** to the amniotic fluid forms of IGFBP-1, **the test can be performed even in the presence of blood and other bodily fluids, infections, and medical products in the extracted specimen.** The high specificity and sensitivity minimize false negative and false positive results, making Actim PROM an optimal solution for diagnosing premature rupture of membranes.

Table 1. Actim PROM has the highest sensitivity, specificity and accuracy in PROM diagnosis (Erdemoglu & Mungan, 2004). When patients with bleeding are also included in the study, the performance of Actim PROM surpasses PAMG-1 tests. (Marcellin et al. 2011)

	Sensitivity	Specificity	Accuracy	
Actim PROM	97 %	97 %	97 %	
pH test	97 %	16 %	56 %	
AFI < 80 mm	94 %	91 %	92 %	

	Sensitivity	Specificity	PPV**	NPV***
Actim PROM	98 %	97 %	98 %	97 %
PAMG-1 test	95 %	95 %	95 %	95 %

Table 2. Clinical evidence of PROM diagnosis with Actim PROM.


	Sensitivity	Specificity	PPV	NPV
Akeran et al., 2005	100 %	92 %	84 %	100 %
Erdemoglu and Mungan, 2004	97 %	97 %	ND	ND
Jain and Morris, 1998	100 %	89 %	76 %	100 %
Rutanen et al., 1996	100 %	95 %	93 %	100 %

*Examples include "Preterm labour and birth", NICE (National Institute for Health and Care Excellence), 2015: www.nice.org.uk/guidance/ng25 (cited on 17.11.2022) and "Premature Birth". Working group established by the Finnish Medical Society Duodecim and the Finnish Gynecological Association. Finnish Medical Society Duodecim, 2018: www.kaypahoito.fi/hoi50089 (cited on 20.09.2022).

**PPV (Positive Predictive Value) = Positive predictive value

***NPV (Negative Predictive Value) = Negative predictive value

Results available at the bedside in just 5 minutes

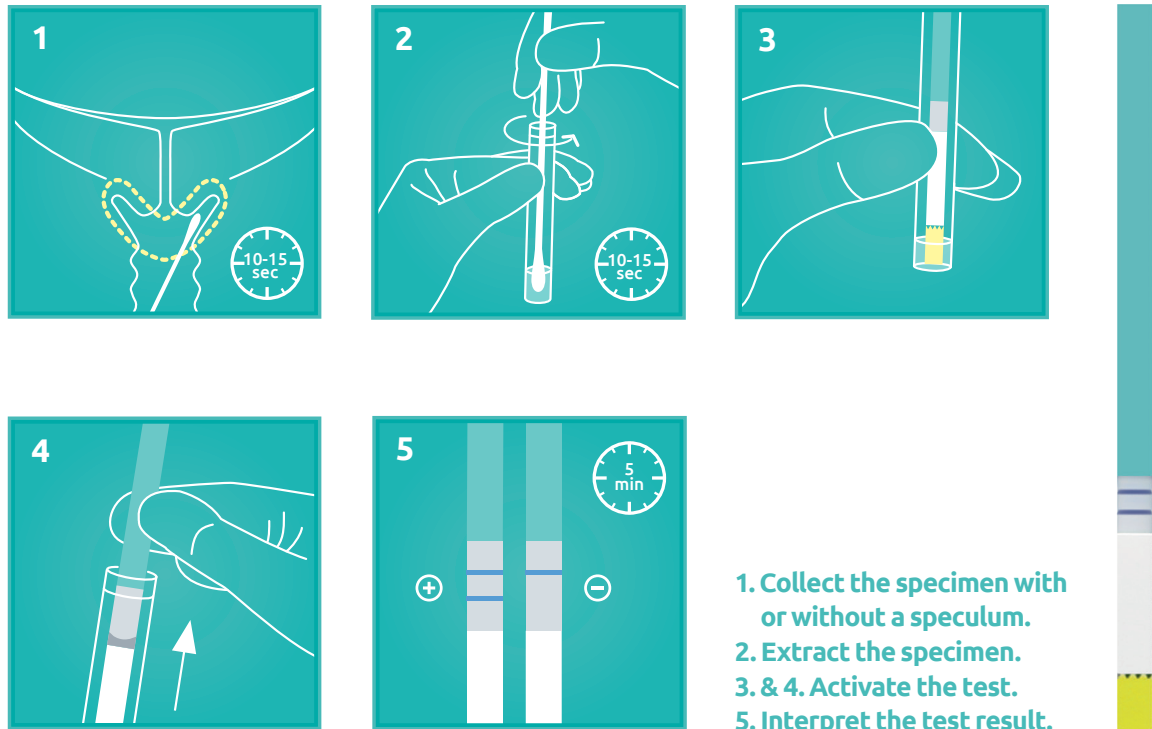

Actim PROM saves time and valuable resources

The diagnosis of PROM is traditionally based on a variety of clinical symptoms and methods. Because symptoms differ among patients, diagnosing this condition can often be challenging and time-consuming.

With Actim PROM, the diagnosis can be made in a timely manner, so treatment can be decided early on to avoid any complications for the pregnant patient and the fetus.

- ✓ Medical attention can be directed to the right patients.
- ✓ Unnecessary use of medication and their side-effects can be avoided.
- ✓ Unnecessary labor inductions can be minimized.
- ✓ Avoidable hospital visits and patient transfers can be reduced.
- ✓ Birth-givers' peace of mind is improved, limiting unnecessary anxiety and worries.

Ordering information

Product	Product code
Actim PROM, 20 tests	30832ETAC
Actim PROM, 10 tests	30831ETAC
Actim PROM, 1 test	30830ETAC
Actim PROM Controls	30800ETAC
Actim 1ngeni Instrument	19101AC
Actim PROM 1ngeni, 10 tests	30831RETAC



The test kit contains all the materials needed - no extra laboratory equipment is required to perform the test.



The test kit can be stored at room temperature.



Actim Oy

– part of Medix Biochemica Group
Klovinpellontie 3, FI-02180 Espoo, Finland
actim@actimtest.com

www.actimtest.com



Selected references

1. Akercan F et al. The value of the insulin-like growth factor binding protein-1 in the cervical-vaginal secretion detected by immunochromatographic dipstick test in the prediction of delivery in women with clinically unconfirmed preterm premature rupture of membranes. *Eur J Obstet Gynecol Reprod Biol* (2005) 121:159-163.
2. Chalurkar U and Andallu R. Detection of Prom using Strip Immuno Assay Test to Detect Insulin Like Growth Factor Binding Protein-1(IGFBP-1) in Amniotic Fluid in Comparison with Fern Test. *Int J Contemporary Medical Research* (2017) 4: 1118-1121.
3. Erdemoglu E and Mungan T. Significance of detecting insulin-like growth factor binding protein-1 in cervicovaginal secretions: Comparison with nitrazine test and amniotic fluid volume assessment. *Acta Obstet Gynecol Scand* (2004) 83: 622-626.
4. Gaucherand P et al. Comparative study of three vaginal markers of the premature rupture of membranes. *Acta Obstet Gynecol Scand* (1997) 76:536-540.
5. Guibourdenche J et al. Rapid detection of insulin-like growth factor-binding protein-1 and foetal fibronectin in cervico-vaginal secretions to diagnose premature membrane rupture. *Ann Clin Biochem* (1999) 36:388-390.
6. Jain K and Morris PG. A clinical study to evaluate the usefulness of the MAST test in diagnosing pre-labour rupture of membranes. *J Obstet Gynaecol* (1998) 18:33-36.
7. Kallioniemi H et al. Usefulness of the insulin-like growth factor binding protein-1 bedside test for ruptured fetal membranes. *Obstet Gynecol Scand* (2014) 93:1282-1289.
8. Kubota T and Takeuchi H. Evaluation of insulin-like growth factor binding protein-1 as a diagnostic tool for rupture of the membranes. *J Obstet Gynecol Res* (1998) 24: 411-417.
9. Marcellin L. et al. Comparison of two bedside tests performed on cervicovaginal fluid to diagnose premature rupture of membranes. *Journal de gynécologie, obstétrique et biologie de la reproduction* (2011) 40:651-656.
10. Palacio M et al. *BMC Pregnancy and Childbirth* (2014), 14:183
11. Rutanen E-M et al. Evaluation of a rapid strip test for insulin-like growth factor binding protein-1 in the diagnosis of ruptured fetal membranes. *Clinica Chimica Acta* (1996) 253: 91-101.
12. Rutanen E-M et al. Measurement of insulin-like growth factor binding protein-1 in cervical/vaginal secretions: comparison with the ROMcheck Membrane immunoassay in the diagnosis of ruptured fetal membranes. *Clinica Chimica Acta* (1993) 214: 73-81.
13. Rutanen E-M. Insulin-like growth factors in obstetrics. *Curr Opin Obstet Gynecol* (2000) 12:163-168.
14. Wathen NC et al. Levels of insulin-like growth factor-binding protein-1 increase rapidly in amniotic fluid from 11 to 16 weeks of pregnancy. *J Endocrinol* (1993) 137:R1-R4.
15. Yang J et al. Vaginal bleeding during pregnancy and preterm birth. *Am J Epidemiol* (2004) 160:118-125.

The full reference list can be found on our website.