



Actim[®]
CRP

Quick CRP
results from
just one drop

In just five minutes Actim[®] CRP rapid test helps estimate if an infection is viral or bacterial. You can get the results anywhere, simply and quickly – with no instrument needed. And all it takes is one small drop of blood.

Our test supports fast decisions on whether your patients require antibiotics or more tests. And the simple sample is especially useful with children and other patients afraid of blood tests.

The Actim CRP rapid test is perfect for outpatient care, or anywhere you need accurate and reliable on the spot results.

About CRP (C-reactive protein)

CRP is an acute-phase protein produced rapidly in response to inflammation, infection, and tissue damage. Level rises very quickly if an infection is present – rising significantly higher in bacterial infections. CRP concentration can also fall very quickly if medical treatment is effective.

Overuse of antibiotics may have devastating consequences. Actim CRP supports your urgent decision on whether the patient requires antibiotics or further tests.


Get CRP results anytime, anywhere

Actim CRP test is extremely easy to use and kit contains everything needed to perform the test. It delivers clear, accurate results in only five minutes. Test requires just 10 µl of whole blood from the fingertip. EDTA, citrate or heparin blood can also be used.

Actim CRP tests can be stored at room temperature, and the tests are conveniently packed in single pouches – so there is no need to take the whole kit to your testing site. Likewise, no laboratory equipment or processing is required, which means zero maintenance or calibrating.

Actim CRP is perfectly suitable for outpatient care or any situation where laboratory services are not available such as schools, nursing homes, workplaces, prisons, and even military camps. The immediate CRP results are also useful in emergency rooms, ambulances, and intensive care units.

Semiquantitative Actim CRP test covers all clinically relevant CRP levels, with three different detection limits: ~10-40 mg/l, 40-80 mg/l and >80 mg/l. With these results, one can provide correct and fast treatment and avoid unnecessarily treating viral infections with antibiotics.

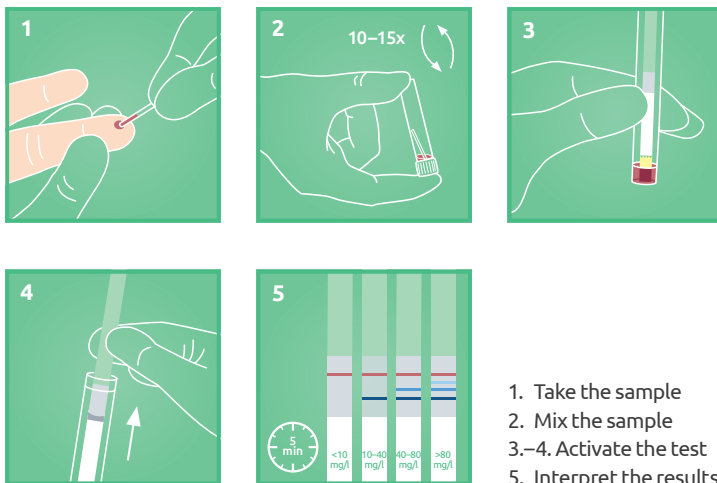

1. Take the sample
2. Mix the sample
- 3.-4. Activate the test
5. Interpret the results

Actim CRP gives results at the bedside in just 5 minutes.

Contact us

Actim Oy

– part of Medix Biochemica Group
Klovinpellontie 3, FI-02180 Espoo, Finland
actim@actimtest.com

www.actimtest.com



Ordering information

Actim CRP 20 test kit

31031ETAC

Product availability and/or labeling may vary from country to country and is subject to varying regulatory requirements. Please contact your local representative for availability. © 10/2021 Actim Oy. All rights reserved. Actim is a registered trademark of Actim Oy.